



EMS Physio Ltd

Wide Neurology Couch

This model is designed for Bobath rehabilitation techniques and the neurological assessment, examination and treatment of both adults and children. The couch features a very wide treatment surface to support both the therapist and patient for these procedures and has an adjustable angle backrest to support the patient in sitting. The couch is available with a choice of hydraulic or electric height operation with push button hand switch control. For very heavy patients the electric model is recommended as it has the capability to raise very heavy loads that cannot be manually achieved on the hydraulic model. A wide range of accessories are available to meet various user requirements.




Standard Features

- UK designed & manufactured
- Lifetime frame & 5 year parts guarantee
- Safe working load electric model 320kg
- Safe working load hydraulic model 250kg
- 16 year average lifespan
- Choice of electric or hydraulic operation
- Breathing hole & plug
- Adjustable angle backrest
- Very strong & rigid frame
- Large treatment surface
- Hardwearing antimicrobial vinyl
- Retractable heavy duty twin wheel base
- Adjustable foot for uneven floors
- Hoist clearance 14cm (5.5")

Overall Dimensions

Length: 189cm (74") approx
 Width: 102cm (40") approx
 Height: 47cm - 97cm approx

Upholstery Dimensions

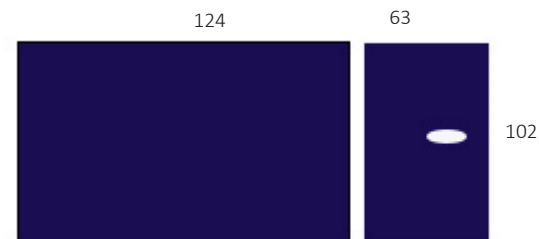

Options & Accessories




Side Security Rails



Paper Roll Holder



Motorised Backrest

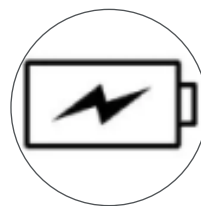

Battery Backup



Footswitch



Large Castor Base

NHS Approved Supplier



t: +44 (0) 1235 772272
 e: sales@emsphysio.co.uk
 www.emsphysio.co.uk

Operational Parameters



Couch at minimum height

EMS Physio Ltd order codes
COU012 Hydraulic with breathing slot
COU008 Electric with breathing slot

Upholstery Vinyl

- Flame retardant to CRIB 5
- Patented anti-microbial technology proven to inhibit the growth of E-Coli and MRSA
- Excellent levels of wear and longevity
- Fully washable and impermeable, to aid infection control
- Suitable for cleaning and disinfecting with dilute solutions of Hypochlorite (10,000ppm)
- Please note; Unless specified Sapphire upholstery will be supplied



Jasmine



Battleship



Cobalt



Cocoa



Grape



Mulled Wine



Pillarbox



Gingersnap



Candy



Applemint



Sapphire



Atlantic Blue



Lupin



Denim



Cool Blue



Rainforest



Wasabi



Citrus Green



Marigold



Almond

For Your Piece of Mind

- Lifetime frame and 5 year parts (accidental damage to upholstery excluded)
- Optional extended warranty to 10 years
- Average life span of 16 years
- Designed and manufactured in-house in the United Kingdom
- CE marked
- Externally audited to ISO9001:2015 by LQRA

Standards

- IP Rating - IPX6
- EE Classification - EU MDR 2017/745 Annex VIII Class 1
- Electrical Standard - IEC 60601-1
- Vinyl flame retardant to CRIB 5 (BS 5852)
- Foam flame retardant to BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service maybe available if required.

EMS PHYSIO Ltd
Downsview Road
Wantage
Oxfordshire OX12 9FE

t: +44 (0) 1235 772272
f: +44 (0) 1235 763518
e: sales@emsphysio.co.uk
www.emsphysio.co.uk